Thermometer Set Exact 70

- | | |
|-------------------------|--|
| INSTRUCTION MANUAL | |
| HANDLEIDING | |
| MANUEL D'UTILISATION | |
| MANUAL DE INSTRUCCIONES | |
| INSTRUKCJA OBSŁUGI | |
| MANUAL DO UTILIZADOR | |

**INTRODUCTION**

Thank you for purchasing the Luvion Digital Thermometer Set Exact 70.

Please retain your original dated sales receipt for your records. For service- and warranty issues you will always be asked for a copy of this receipt. For warranty issues you can contact your local Luvion dealer or distributor. More information can be found on www.luvion.com.

The user manual provides you with all information for a safe and easy use of this thermometer set. Before using the thermometer set we ask you to thoroughly read this user manual. We ask for your special attention for the safety instructions.

CONTENTS OF PACKAGE

- Digital flex-tip thermometer
- Digital soother thermometer
- Digital bath and room thermometer

SAFETY INSTRUCTIONS

- Only use the flex-tip and soother thermometers to measure human body temperature.
- Before each use check the nipple of the soother for tears or cracks and pull it to ensure that it is firmly attached to the shield. Do not use it if a defect is found or if it is loose.
- Ensure an adult supervises the use of the thermometers.
- Store the thermometers out of reach of children; they are not toys.
- Temperature measurement is not a substitute for consulting your physician.

CARING FOR THE THERMOMETERS

- Do not store the thermometers in direct sunlight or at a high temperature.

- Do not disassemble the thermometers.
- Never use boiling water or a chemical solvent to clean the thermometers.
- Do not wash in dishwasher.
- Do not put in microwave or microwave steriliser.

BEFORE FIRST USE

Clean the flex-tip thermometer and the soother thermometer before first use.

THE DIGITAL FLEX-TIP THERMOMETER

The digital flex-tip thermometer can be used to measure your baby's temperature anally (rectum) or axially (armpit).

USING THE FLEX-TIP THERMOMETER

1. Remove the protection cap from the thermometer.
2. Press the power button to switch the thermometer on. The display will read "188.8E," confirming that the thermometer is working.
3. Choose the place of the temperature measurement. Anal measurements are most accurate and are suitable for both infants and children. Axial measurements are not recommended for accurate measurements.
 - When used anally (rectum), the temperature sensor needs to be inserted for about 1 to 2 cm. Vasaline can be used for lubrication.
Measurement time 30 seconds
Normal temperature about 37°C-98.6°F

- When used axially (armpit), put the temperature sensor in the armpit. Press and fix it with the arm.
Measurement time 3~5 minutes
Normal temperature about 36.7°C-98.1°F

⚠ Note: The body temperature is different for all humans and changes during the day. You are therefore advised to measure the temperature at a specific time in the morning.

- Start the measurement when the screen displays "Lo °C" or "Lo °F". During measurement the symbol "°C" or "°F" will flash while the temperature increases. When a maximum value is reached, the symbol "°C" or "°F" will cease to flash. A buzzing sound that lasts about 10 seconds indicates that the measurement is completed.

10 beeps in 10 seconds means that the measured temperature is normal (below 37.8°C) 3 sets of short beeps in 10 seconds means a possible fever (measured temperature above 37.8°C).

- When the display reads "Lo °C" or "Lo °F" then the measured temperature is below 32.0°C/ 90 °F
- When the display reads "Hi °C" or "Hi °F" then the measured temperature is above 43.0°C/ 109.9 °F

- Press the power button to turn the thermometer off.

⚠ Note: The thermometer will automatically switch off after 10 minutes. Switch the thermometer off directly after use to prolong battery life.

RECALL THE LAST TEMPERATURE

Press the power button while turning on the flex-tip thermometer. The display will now read the last measured temperature. This is indicated by the letter "M" displayed under the degree sign. After 3 seconds the display will read 37.0°C and is ready for a new measurement.

CLEANING

The thermometer is waterproof and can be cleaned with water, sterilization lotion or a soft cloth.

REPLACING THE BATTERY

When the display becomes indistinct or when the “” symbol appears it is time to replace the battery.

To replace the battery:

1. Open the battery cover, the blue cover at the top of the thermometer.
2. Remove the old battery.
3. Place a new battery and reattach the battery cover.

SPECIFICATIONS

Measuring Range	32.0°C ~ 42.9°C (89.6°F ~ 109.2°F)
Measuring Time	About 3 ~ 5 minute
Accuracy	± 0.1°C (±0.2°C) for 35.0 ~ 42.0°C (95°F ~ 107.6°F) & ±0.2°C (±0.4°F) for remaining range at an ambient of 22°C (71.6°F)
Display	LCD (Unit of display : 0.1°C / 0.1°F)
Battery	1.5VDC (LR41) alkaline
Power Consumption	0.15 mW
Dimensions	129mm(L) x 21mm(W) x 12.5mm(H)
Battery Life	About 250 hours or 1000 times.
Weight	Approx. 12 g. (with battery)

THE SOOTHER THERMOMETER

The most convenient way to measure your baby's temperature is with a soother thermometer.

Note: Oral (mouth) measurements are often about 0.5 °C or 1.0 °F lower than body temperature readings (anal measurement).

USING THE SOOTHER THERMOMETER

1. Remove the protection cap from the soother.
2. Press the power button. The display will read "188.8E," confirming that the thermometer is working.
3. Place the soother in the baby's mouth.

The measurement starts when the screen displays "Lo °C" or "Lo °F". During measurement the symbol "°C" or "°F" will flash while the temperature increases. When a maximum value is reached (after approx. 3 minutes), the symbol "°C" or "°F" will cease to flash. A buzzing sound indicates that the measurement is finished.

Note: It is recommended to leave the soother thermometer in the baby's mouth for 6-8 minutes to allow the soother to warm up and to get a more accurate measurement.

- When the display reads "Lo °C" or "Lo °F" then the measured temperature is below 32.0°C/ 90 °F
- When the display reads "Hi °C" or "Hi °F" then the measured temperature is above 43.0°C/ 109.9 °F

Note: Drinking or eating prior to use may affect the temperature reading.
Wait half an hour before taking a measurement.

4. Press the power button to switch off the soother thermometer.

Note: The thermometer will automatically switch off after 10 minutes.
Switch the thermometer off directly after use to prolong battery life.

RECALL THE LAST TEMPERATURE

To recall the last temperature:

Press and hold the power button while turning on the soother thermometer

for at least 2 seconds. The display will now read the last measured temperature. This is indicated by the letter "M" displayed under the degree sign.

CLEANING

Wash the soother thermometer after use with luke warm water and a mild soap.

BATTERY

When the battery is almost exhausted, the "■" symbol will appear on the display. This means that it is time to discard the unit.

 **Note:** Do not use the soother thermometer as a regular soother after the battery is exhausted.

SPECIFICATIONS

Measuring Range	32.0°C ~ 43.9°C (89.6°F ~ 109.9°F)
Accuracy	± 0.1°C (±0.2°C) for 35.0 ~ 38.0°C (95°F ~ 100.4°F) ambient temperature of 22°C (71.6°F)
Display	LCD (Unit of display : 0.1°C / 0.1°F)
Power Consumption	0.15 mW
Dimensions	40mm(L) x 43mm(W) x 53mm(H)
Battery Life	More than 200 hours of continuous use. (approx. 2 years if used 10 minutes a day)
Weight	Approx. 14 g. (with battery)

THE BATH AND ROOM THERMOMETER

Use the Bath and Room Thermometer to conveniently measure the temperature of the baby room or the bathwater.

Note: Your baby will feel most comfortable in bath if the water temperature is between 36.5°C and 38°C.

! A temperature of 39°C and above is too hot and your baby could get burned! At a room temperature of about 18°C babies feel most comfortable when sleeping.

USING THE BATH AND ROOM THERMOMETER

Switching on: Switch on the bath thermometer by firmly hitting the bottom side with a flat hand.

Just let the thermometer float in the bathwater or let it stand on a shelf in the baby's room to start measuring.

Note: Mix the water thoroughly before measuring the water temperature

! or put the thermometer in the tub while filling it. This will ensure a fast and accurate reading of the water temperature.


- The display will read "COLD" when the measured temperature is below 33°C. A comfortable bathwater temperature is between 36.5°C and 38°C.
- Warm water alert: the red warning LED will start to flash when the measured temperature is 39 °C or above. The display will then flash "Hot."

REPLACING THE BATTERY

After about 2 years the batteries may go flat. You can then replace them.

To replace the batteries:

1. Push the measurement unit out of the rubber housing
2. Remove the battery cover located at the back of the measurement unit by using a coin or similar
3. Carefully remove the empty batteries with a pair of tweezers

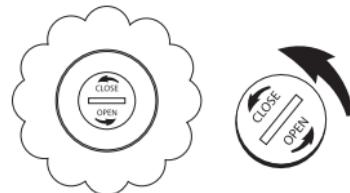SPECIFICATIONS

Measuring Range	0°C ~ 50°C (32°F ~ 122°F)
Accuracy	± 1°C (2°F)
Display	LCD (Unit of display : 0.1°C / 0.1°F)
Battery	1x CR2032
Dimensions	40mm(L) x 43mm(W) x 53mm(H)
Battery Life	about 1 year (based on 1 hour of use a day)

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